



TRIANGO 80 LED PROCEDURE LIGHT



PRODUCT INFORMATION: TRIANGO 80

TRIANGO 80 C

LED technology has revolutionized the art of equipment design for the medical sector.

TRIANGO 80 brings this innovative technology to the procedure & exam area. With its advanced thermal management system, integrated dimming controls, brilliant light output, and compact housing, TRIANGO 80 meets the stringent requirements of both hospital designers and practitioners.




**Uniform illumination
for work area**

High intensity light
80,000 lux / 1 meter

Glare-free
To combat eye fatigue

Easy to clean
The removable handle can be sterilized
and the closed construction of the
fixture allows for easy cleaning

Unparalleled color rendering
Color Rendering Index (CRI) > 95

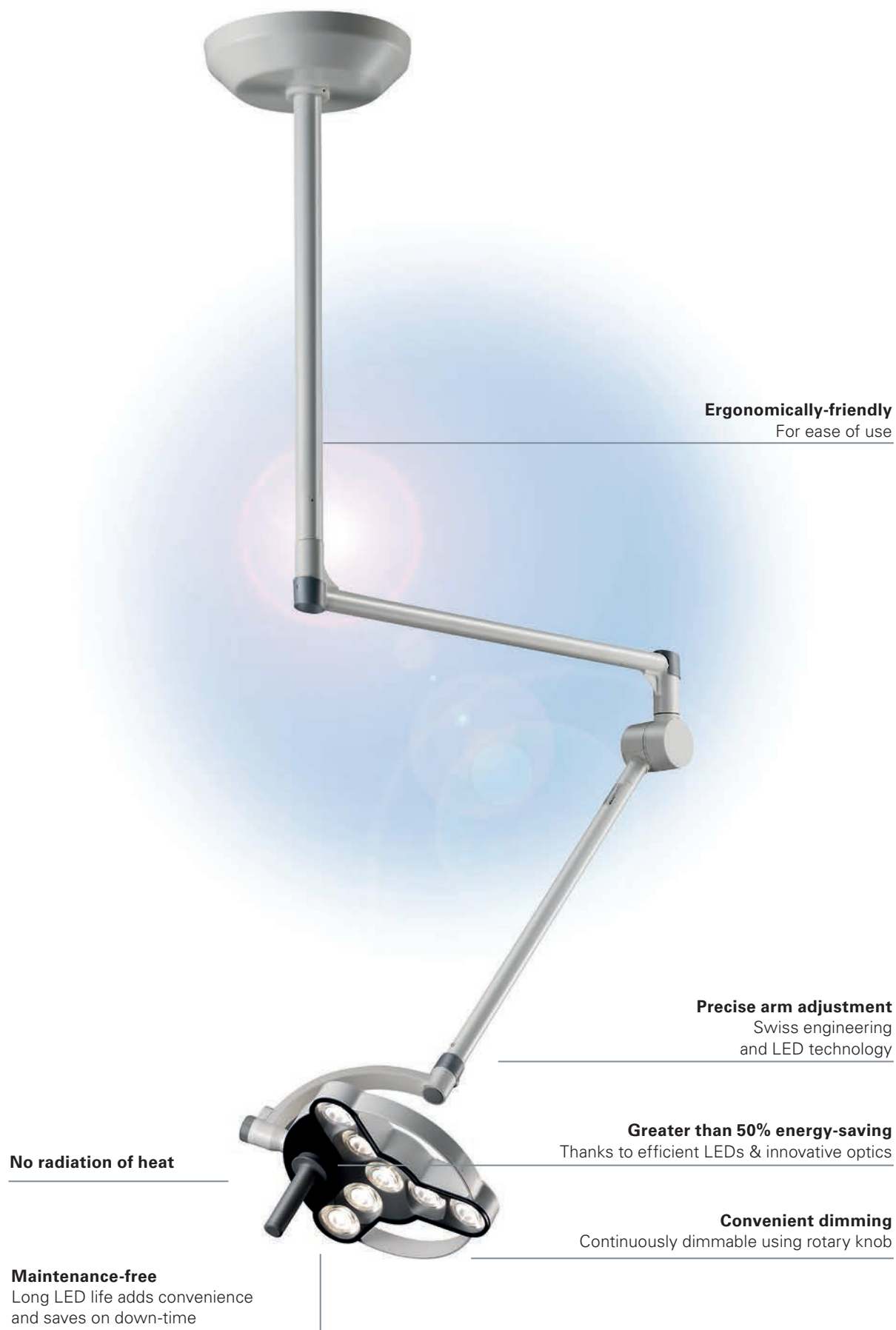



PRODUCT ADVANTAGES

TRIANGO 80



MOBILITY AND FUNCTION

LED technology gives high performance lighting for procedure areas. TRIANGO 80's intricate design and robust materials also allow for a broad range of movement and easy positioning of the fixture.

AVAILABLE IN MULTIPLE VERSIONS

Ceiling Mount:


The ceiling mount can be a space saver, especially in rooms where maximum freedom of movement is required. Ceiling mount model is available in single-head version.

Floor Stand:

Taking the TRIANGO 80 mobile with the rolling floor stand offers flexibility for various areas of the room or even sharing among several rooms.



Wall Mount:

The wall mounted version is also available.

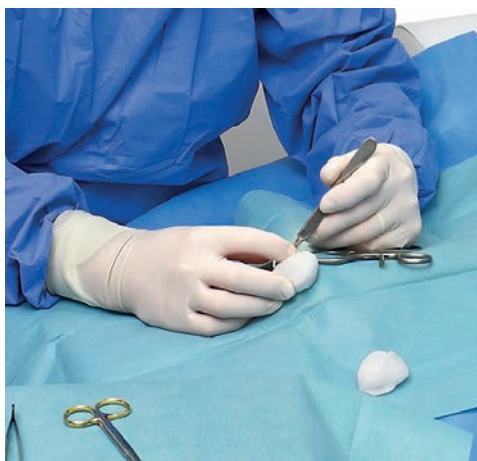

EVEN ILLUMINATION OF THE AREA

The unique optic developed by Derungs ensures an optimum light field. The illuminated area meets the highest technical lighting requirements for procedure lights. The reflection-, shadow- and glare-free lighting generates an ideal field of vision specific to this type of close work.

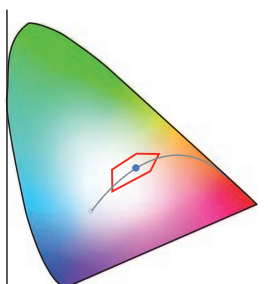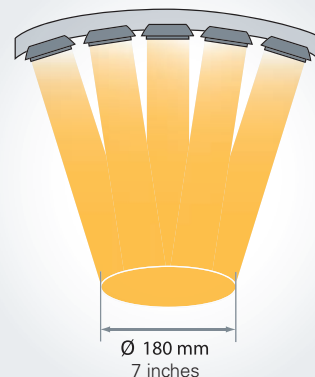


- Black Body Curve
- EN 60601-2-41
- ● ● optimal light color temperature

COLOR RENDERING

TRIANGO 80 offers the optimum in color rendering with $R_a > 95$ and $R_9 > 90$. This allows for clear recognition of any variations or irregularities in tissue while performing procedures.



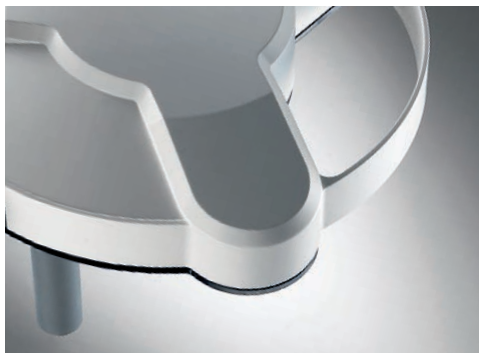

MAINTENANCE-FREE

LED technology allows for long periods of maintenance-free usage. No service for bulb replacement necessary for the life of the fixture.

- Long lamp life of 50,000 hours

CLEANING AND STERILIZATION

Its closed construction and smooth surface allows for quick, thorough cleaning. In addition, any penetration of dust, dirt or liquids is prevented. The removable handle is easy to sterilize.

- Removable and sterilizable handle
- Closed construction protects mechanics



HIGH-QUALITY COMPONENTS

High-quality components guarantee a strong, user-friendly luminaire. Derungs uses only top of the line materials to garner maximum benefits from its product.

Luminaire body:

Polycarbonate & acrylonitrile-butadiene-styrene copolymer (PC & ABS)

Lens Cover:

Polymethyl methacrylate (PMMA)

An ABR-GL protective coating makes the lens cover scratch proof.

IMPROVED ENERGY EFFICIENCY

Notable energy savings is gained thanks to the most up-to-date LED technology: high illumination using only 30 watts of power.

If lower light levels are necessary, power consumption will be even less.



ILLUMINATION / LIGHT INTENSITY REGULATION

Easy operation: One control switches the light on/off and dims with five light levels, providing the practitioner with the quick, convenient control necessary when many different procedures are being done during the course of the day.

GLARE-FREE

The light beam is focused and glare-free whatever the angle of sight. This avoids bothersome reflections on the surface.

- Continuous dimming
- 80,000 lux/39.5 in (1m)
- Illuminated area 7 in (180mm)
- CRI > 95
- Glare-, shadow- and reflection-free



NO RADIATION OF HEAT

An ideal working environment for a physician during a procedure includes clear, glare-free light with almost no infrared or UV radiation. Not only does this add to the comfort of the physician, staff and patient, but also prevents tissue from drying out or being affected by the heat, a genuine concern to the provider.

GOOD WORKING CONDITIONS

Out of necessity, any procedure light is often positioned in close proximity to patient and staff. Optics are used to create an environment free of shadows, reflections and glare.

Modern LEDs and innovative design and technology make TRIANGO one of the most effective tools used in this arena.

ORDERING INFORMATION

Ceiling Mount

Part No	Mounting Type
D16138000	Single Ceiling

Ceiling Mount



Wall Mount

Part No	Mounting Type
D16139000	Wall

Wall Mount



Floor Stand

Part No	Mounting Type
D16135000	Floor Stand

Floor Stand



HEADQUARTERS GERMANY
Herbert Waldmann GmbH & Co. KG
Peter-Henlein-Straße 5
78056 VILLINGEN-SCHWENNINGEN
GERMANY
Telephone +49 7720 601-0
Telephone +49 7720 601-100 (Sales)
Fax +49 7720 601-290
www.waldmann.com
sales.germany@waldmann.com

FRANCE
Waldmann Eclairage S.A.S.
Zone Industrielle
Rue de l'Embranchement
67116 REICHSTETT
FRANCE
Telephone +33 3 88 20 95 88
Fax +33 3 88 20 95 68
www.waldmann.com
info-fr@waldmann.com

ITALY
Waldmann Illuminotecnica S.r.l.
Via della Pace, 18 A
20098 SAN GIULIANO MILANESE (MI)
ITALY
Telephone +39 02 98 24 90 24
Fax +39 02 98 24 63 78
www.waldmann.com
info-it@waldmann.com

THE NETHERLANDS
Waldmann B.V.
Lingewei 19
4004 LK TIEL
THE NETHERLANDS
Telephone +31 344 631019
Fax +31 344 627856
www.waldmann.com
info-nl@waldmann.com

AUSTRIA
Waldmann Lichttechnik Ges.m.b.H.
Gewerbepark Wagram 7
4061 PASCHING/LINZ
AUSTRIA
Telephone +43 7229 67400
Fax +43 7229 67444
www.waldmann.com
info-at@waldmann.com

SCANDINAVIA
Waldmann Ljusteknik AB
Skebokvarnsvägen 370
124 50 BANDHAGEN
SWEDEN
Telephone +46 8 990 350
Fax +46 8 991 609
www.waldmann.com
info-se@waldmann.com

SWITZERLAND
Waldmann Lichttechnik GmbH
Benkenstrasse 57
5024 KÜTTIGEN
SWITZERLAND
Telephone +41 62 839 1212
Fax +41 62 839 1299
www.waldmann.com
info-ch@waldmann.com

UNITED KINGDOM
Waldmann Lighting Ltd.
10 Millfield House
Croxley Park
WATFORD WD18 8YX
UNITED KINGDOM
Telephone +44 1923 800030
Fax +44 1923 800016
www.waldmann.com
info-uk@waldmann.com

USA
Waldmann Lighting Company
9 W. Century Drive
WHEELING, ILLINOIS 60090
USA
Telephone +1 847 520 1060
Fax +1 847 520 1730
www.waldmannlighting.com
waldmann@waldmannlighting.com

CHINA
Waldmann Lighting (Shanghai) Co., Ltd.
Part A11a, No. Five Normative Workshop
199 Changjian Road, Baoshan
SHANGHAI, P.R.C. 200949
CHINA
Telephone +86 21 5169 1799
Fax +86 21 3385 0032
www.waldmann.com.cn
info-cn@waldmann.com

SINGAPORE
Waldmann Lighting Singapore Pte. Ltd.
77A Neil Road
SINGAPORE 088903
SINGAPORE
Telephone +65 6275 8300
Fax +65 6275 8377
www.waldmann.com
sales-sg@waldmann.com