

Varios Combi Pro2

Combo Piezo Scaler & Air Polisher



**ULTRASONIC
SCALER**

**POWDER
THERAPY**

Combining subgingival & supragingival powder therapy
with ultrasonic scaling for periodontal treatment

Constant Availability



Independent Unit & Chambers

The powder chamber, hose, and handpieces are separate from the main unit. In the unlikely event of powder clogging, the parts can easily be detached for cleaning or service.

Patient Comfort



Water Heating System

Injected water is warmed to approximately the same temperature as the oral cavity to provide ultimate comfort, even those with sensitivity.

Easy Care



Auto-Cleaning

The 4 auto-cleaning modes help prevent clogging and reduce water accumulation.



CLM Mode

Easily reprocess the equipment using the CLM Mode and cleaning solution for a hands-free cleaning.

Operability



Wireless Foot Control

The multi-functional wireless foot control allows for optimal positioning.

Varios Combi Pro2

Basic set



OPTIC	MODEL	ORDER CODE
●	Varios Combi Pro2 Basic Set	Y1501049

Contents:

- Varios Combi Pro2 Control Unit
- VA2-LUX-HP Ti
- VA2-LUX-SC Ti
- Prophy-Chamber2
- Prophy-Hose2
- Powder Handpiece (Prophy)
- Wireless Foot Control
- Hanger
- Air Hose
- Water Hose
- Power cord
- Irrigation Bottle Set
- CLM Bottle Set
- Ultrasonic Scaler Tips (G8, G16, P20)
- FLASH pearl powder (300g)

Specifications:

- Dimensions: W 249 x D 290 x H 115 mm
- Weight: 3.5 kg
- Power source: AC 120 V 50/60 Hz
- Drive frequency (ultrasonic): 28 - 32 kHz
- Maximum output (ultrasonic): 11 W



Scan to
download
Tip Guide

Add-On



MODEL	ORDER CODE
Perio Kit2	Y1500470

Contents:

- Perio-Chamber2
- Perio-Hose2
- Powder Handpiece (Perio)
- Perio Nozzle Tip 40pc
- Perio Mate Nozzle tip remover
- Perio Mate Powder (160g)

FLASH pearl

Polishing powder for supragingival use, with spherical particles that are gentle on the tooth surface.
Intended use: removal of stubborn stains




54µm
Calcium
Carbonate

MODEL	ORDER CODE
FLASH pearl 300g bottle x 4	Y900698

Perio Mate Powder

Polishing powder for gingival margin & subgingival use.
Intended use: removal of biofilm and minor staining

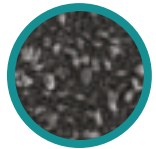

25µm
Glycine

MODEL	ORDER CODE
Perio Mate Powder 160g bottle x 2	Y900938

Powder Therapy Combinations & Uses

MODE	KIT	HANDPIECE	NOZZLE	USE	POWDER
PROPHY 	Prophy Kit2 				
	Perio Kit2 				
PERIO 	Perio Kit2 				
	Perio Kit2 				